



VirtuMED

OKAY II™


Y-Connector Set



Y-Connector

 **NIPRO**

Features

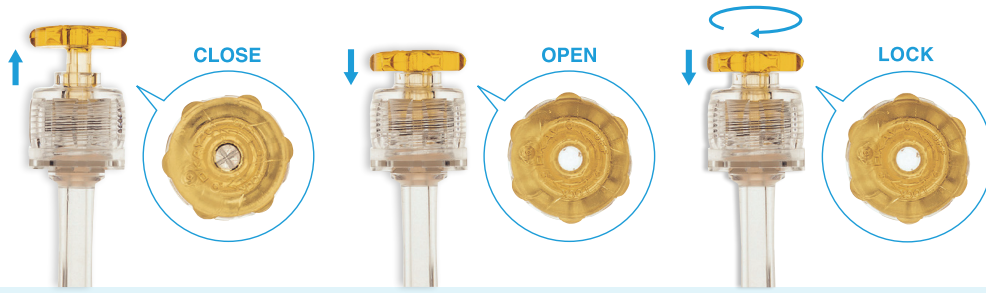

Hemostatic valve

The hemostatic valve can be opened by pushing the opener, and closed by releasing it.

Once the hemostatic valve has been opened, a clockwise quarter turn locks the valve in the open position.



Fixed valve

Valve to fix devices and prevent loss of pressure.

Large lumen

Nominal ID: 10F, 3.33mm

Allows for use in combination with various products for a variety of PCI strategies.

Superior design

Facilitates single-handed operation.


Ease-of-use

Hemostatic valve features a locking mechanism.

Minimizes blood loss

Two valves (hemostatic and fixed) greatly reduce blood loss during the procedure.

Accessories

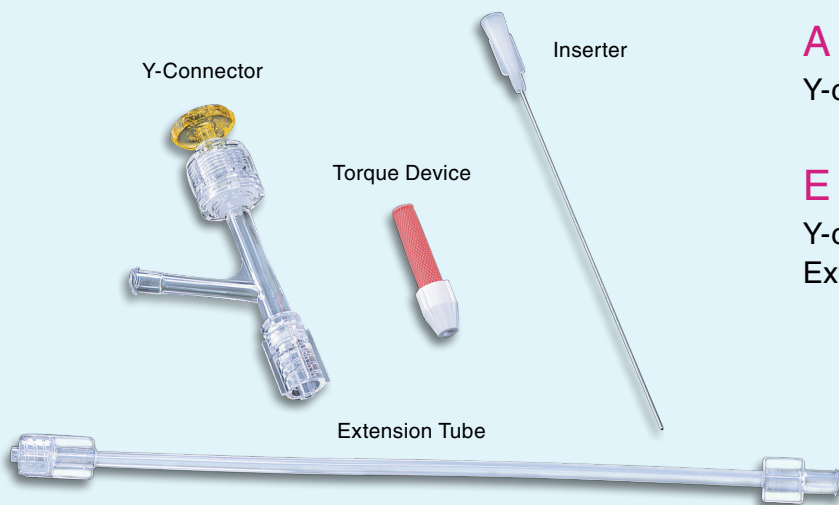

A set (YOK0A)

Y-connector, Inserter, Torque Device

E set (YOK0E)

Y-connector, Inserter, Torque Device, Extension tube